


Is your SpiroCheck Duo ready for a tune-up?

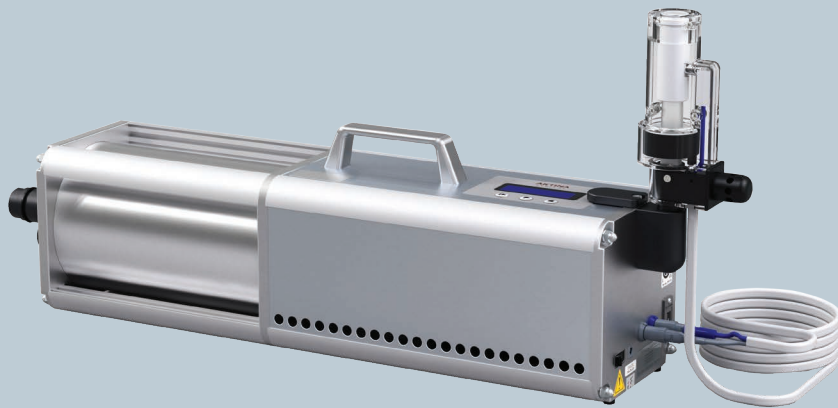

Did you know the SpiroCheck Duo ABC QA tool is shipped from the Aktina Medical factory with a 3-year guaranteed volume?

Are you due for a calibration?

For a quick answer have a look at the calibration sticker on the front of the unit (see Figure 1). The calibration timeframe is clearly indicated. If you need a calibration service, you will also benefit from the Aktina Medical preventive maintenance which is included.

Aktina Medical		Calibration Summary Tag	
Model: 12-202 SpiroCheck Duo			
APR 17, 2019			
Cal Date	2050	2.998L	
APR 17, 2022	Serial	Calibrated Volume	
Cal Due			
(888) 433-3380 info@aktina.com			

Figure 1



The calibration service includes a high accuracy measurement of the SpiroCheck Duo unit for displaced air volume. Technicians then update the internal software to incorporate the true measured volume which is required for pass/fail determination by the internal software. We will also perform a comprehensive preventive maintenance check on the full system and make corrections as necessary to pneumatic tubing, fittings, and gaskets. Note: this does not include replacement of any internal hardware.

This comprehensive work is usually performed within a one-week turnaround time.

Please contact Aktina Medical at info@aktina.com or 888.433.3380 for more information.

